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BABYWAAGE 100KG

USER MANUAL

BABY SCALE 100KG

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PÈSE-BÉBÉ 100KG

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BILANCIA NEONATO 100KG

MANUAL DE INSTRUCCIONES

BALANZA PARA BEBÉ 100KG

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SBS-BS-100C

▪ Deutsch	3
▪ English	5
▪ Polski	7
▪ Česky	9
▪ Français	11
▪ Italiano	13
▪ Español	15

PRODUKTNAME	BABYWAAGE 100KG
PRODUCT NAME	BABY SCALE 100KG
NAZWA PRODUKTU	WAGA DZIECIĘCA 100KG
NÁZEV VÝROBKU	DĚTSKÁ VÁHA 100KG
NOM DU PRODUIT	PÈSE-BÉBÉ 100KG
NOME DEL PRODOTTO	BILANCIA NEONATO 100KG
NOMBRE DEL PRODUCTO	BALANZA PARA BEBÉ 100KG
MODELL	
PRODUCT MODEL	
MODEL PRODUKTU	
MODEL VÝROBKU	SBS-BS-100C
MODÈLE	
MODELLO	
MODELO	
HERSTELLER	
MANUFACTURER	
PRODUCENT	
VÝROBCE	EXPONDO POLSKA SP. Z O.O. SP. K.
FABRICANT	
PRODUTTORE	
FABRICANTE	
ANSCHRIFT DES HERSTELLERS	
MANUFACTURER ADDRESS	
ADRES PRODUCENTA	
ADRESA VÝROBCE	UL. NOWY KISIELIN-INNOWACYJNA 7, 66-002 ZIELONA GÓRA POLAND, EU
ADRESSE DU FABRICANT	
INDIRIZZO DEL PRODUTTORE	
DIRECCIÓN DEL FABRICANTE	

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Parameter	Werte
Produktname	BABYWAAGE 100KG
Modell	SBS-BS-100C
Abmessungen [mm]	355x308x30
Gewicht [kg]	1,8
Mindestlast [g]	50
Maximallast [kg]	100
Genauigkeit [g]	<40 kg 10 >40kg 100
Batterieparameter	1,5V 3xAAA
Automatische Abschaltzeit [s]	300

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1. GENERAL DESCRIPTION

The user manual is designed to assist in the safe and trouble-free use of the device. The product is designed and manufactured in accordance with strict technical guidelines, using state-of-the-art technologies and components. Additionally, it is produced in compliance with the most stringent quality standards.

DO NOT USE THE DEVICE UNLESS YOU HAVE THOROUGHLY READ AND UNDERSTOOD THIS USER MANUAL.

To increase the product life of the device and to ensure a trouble-free operation, use it in accordance with this user manual and regularly perform the maintenance tasks. The technical data and specifications in this user manual are up-to-date. The manufacturer reserves the right to make changes associated with quality improvement.

LEGEND

- The product satisfies the relevant safety standards.
- Read instructions before use.
- The product must be recycled.
- WARNING!** or **CAUTION!** or **REMEMBER!** Applicable to the given situation (general warning sign).
- Only use indoors.

PLEASE NOTE! Drawings in this manual are for illustration purposes only and in some details may differ from the actual machine.

The original operation manual is in German. Other language versions are translations from German.

2. USAGE SAFETY

ATTENTION! Read all safety warnings and all instructions. Failure to follow the warnings and instructions may result in an electric shock, fire and/or serious injury or even death.

The terms "device" or "product" are used in the warnings and instructions to refer to BABY SCALE 100KG.

2.1. SAFETY IN THE WORKPLACE

- a) Make sure the workplace is clean and well lit. A messy or poorly lit workplace may lead to accidents. Try to think ahead, observe what is going on and use common sense when working with the device.
- b) If you are unsure about whether the product is operating correctly or if you find damage, please contact the manufacturer's service centre.
- c) Only the manufacturer's service centre may make repairs to the product. Do not attempt to make repairs yourself!
- d) Please keep this manual available for future reference. If this device is passed on to a third party, the manual must be passed on with it.
- e) Keep packaging elements and small assembly parts in a place not available to children.

2.2. PERSONAL SAFETY

- a) The device is not designed to be handled by persons (including children) with limited mental and sensory functions or persons lacking relevant experience and/or knowledge unless they are supervised by a person responsible for their safety or they have received instruction on how to operate the device.
- b) The device can be handled only by physically fit persons who are capable of handling it, properly trained, familiar with this manual and trained within the scope of occupational health and safety.
- c) The device is not a toy. Children must be supervised to ensure that they do not play with the device.

2.3. SAFE DEVICE USE


- a) When not in use, store in a safe place, away from children and people not familiar with the device who have not read the user manual.
- b) Keep the device out of the reach of children.
- c) Device repair or maintenance should be carried out by qualified persons, only using original spare parts. This will ensure safe use.
- d) To ensure the operational integrity of the device, do not remove factory-fitted guards and do not loosen any screws.
- e) Clean the device regularly to prevent stubborn grime from accumulating.
- f) The device is not a toy. Cleaning and maintenance may not be carried out by children without supervision by an adult person.
- g) It is forbidden to interfere with the structure of the device in order to change its parameters or construction.
- h) Keep the device away from sources of fire and heat.
- i) Do not overload the device.

ATTENTION! Despite the safe design of the device and its protective features, and despite the use of additional elements protecting the operator, there is still a slight risk of accident or injury when using the device. Stay alert and use common sense when using the device.

3. USE GUIDELINES

The device is used for precise weight measurement. **The user is liable for any damage resulting from improper use of the device.**

3.1. DEVICE DESCRIPTION

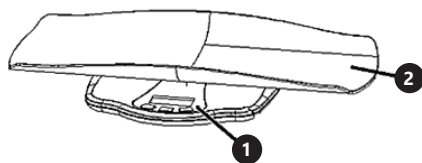1. Control panel with display
2. Scale pan

3.2. DEVICE USE

3.2.1. ASSEMBLING THE DEVICE

⚠ ATTENTION! This product's exploded view can be found on the last pages of the operating instructions (pp. 17).

3.2.2. CONTROL PANEL

- 1) T button - it is used to turn the device on or off (when the pan is empty) and for taring (when there is a child or weights on the pan)
- 2) UNIT button: used to change the measurement unit
- 3) HOLD button : used to lock the reading value (the value will not disappear from the display even though a product is removed from the pan)

3.2.3. WEIGHTING

1. Place the scale on an even, stable, dry and hard surface, away from any hot surfaces and out of the reach of children.
2. Insert the batteries into the device
3. Press the T button.
4. Wait until the display shows „0“.
5. Put a child on a scale.
6. Read the value after a few seconds when the measurement is stabilized (press UNIT to change the unit).

CAUTION: The display may show the indicator:

- „LO“ - informs about a low battery level,
- „O_Ld“ - informs about overload,
- UnSt- means instability (it may appear if a child is moving during weighing, then weighing should be repeated).

CAUTION: If the device is not used, it will turn off automatically after 5 minutes.

3.2.4. CALIBRATION

1. Place the scale on an even, stable and hard surface.
2. Press the button: 3 times briefly and 1 time long until the display shows the internal code.
3. Press the button briefly again.
4. The display will show „25“, put the weight of 25 kg on the pan.
5. Then the display will show „50“, place a 50 kg weight on the pan.
6. Then the display will show „75“, place a 75 kg weight on the pan.
7. Then the display will show „100“, place a 100 kg weight on the pan.
8. Place the appropriate weight (according to the value selected) on the pan. When the (last) reading is stabilized, the display will show „PASS“.

CAUTION: If the weight becomes unstable or the calibration is interrupted, the display will show "Err" and will turn off automatically.

9. Turn the scale off.

3.3. CLEANING AND MAINTENANCE

- a) Always unplug the device before cleaning it.
- b) Store the unit in a dry, cool place, free from moisture and direct exposure to sunlight.
- c) Do not spray the device with a water jet or submerge it in water.
- d) The device must be regularly inspected to check its technical efficiency and spot any damage.
- e) Use a soft cloth for cleaning.
- f) Do not leave the battery in the device if it will not be used for a longer period of time.
- g) Do not use sharp and/or metal objects for cleaning (e.g. a wire brush or a metal spatula) because they may damage the surface material of the appliance.
- h) Replace the batteries when the display shows „Lo“.

SAFE REMOVAL OF BATTERIES AND RECHARGEABLE BATTERIES

1.5V AAA batteries are installed in the devices. Remove used batteries from the device using the same procedure by which you installed them. Recycle batteries with the appropriate organisation or company.

DISPOSING OF USED DEVICES

Do not dispose of this device in municipal waste systems. Hand it over to an electric and electrical device recycling and collection point. Check the symbol on the product, instruction manual and packaging. The plastics used to construct the device can be recycled in accordance with their markings. By choosing to recycle you are making a significant contribution to the protection of our environment. Contact local authorities for information on your local recycling facility.

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