USER MANUAL

expondo.com

BEDIENUNGSANLEITUNG
INSTRUKCJA OBSŁUGI
NÁVOD K POUŽITÍ
MANUEL D'UTILISATION
ISTRUZIONI PER L'USO
MANUAL DE INSTRUCCIONES
HASZNÁLATI ÚTMUTATÓ
BRUGSANVISNING

PHY_SMC_01

MASSAGE CHAIR

DE	Produktname	MASSAGESESEL
EN	Product name	MASSAGE CHAIR
PL	Nazwa produktu	FOTEL MASUJĄCY
CZ	Název výrobku	MASÁŽNÍ KŘESLO
FR	Nom du produit	FAUTEUIL DE MASSAGE
IT	Nome del prodotto	POLTRONA MASSAGGIANTE
ES	Nombre del producto	SILLÓN DE MASAJE
HU	Termék neve	MASSÁZSFOTEL
DA	Produktnavn	MASSAGESTOL
FI	Tuotteen nimi	HIERONTATUOLI
NL	Productnaam	MASSAGESTOEL
NO	Produktnavn	MASSASJESTOL
SE	Produktnamn	MASSAGESTOL
PT	Nome do produto	CADEIRA DE MASSAGEM
SK	Názov produktu	MASÁŽNE KRESLO
DE	Modell	PHY_SMC_01
EN	Product model	
PL	Model produktu	
CZ	Model výrobku	
FR	Modèle	
IT	Modello	
ES	Modelo	
HU	Modell	
DA	Model	
FI	Tuotteen malli	
NL	Productmodel	
NO	Produktmodell	
SE	Produktmodell	
PT	Modelo do produto	
SK	Model	
DE	Hersteller	expondo Polska sp. z o.o. sp. k.
EN	Manufacturer	
PL	Producent	
CZ	Výrobce	
FR	Fabricant	
IT	Produttore	
ES	Fabricante	
HU	Termelő	
DA	Producent	
FI	Valmistaja	
NL	Producent	
NO	Produsent	
SE	Tillverkare	
PT	Fabricante	
SK	Výrobca	
DE	Anschrift des Herstellers	ul. Nowy Kisielin – Innowacyjna 7, 66-002 Zielona Góra Poland, EU
EN	Manufacturer Address	
PL	Adres producenta	
CZ	Adresa výrobce	
FR	Adresse du fabricant	
IT	Indirizzo del produttore	
ES	Dirección del fabricante	
HU	A gyártó címe	
DA	Producentens adresse	
FI	Valmistajan osoite	
NL	Adres producent	
NO	Produsentens adresse	
SE	Tillverkarens adress	
PT	Endereço do fabricante	
SK	Adresa výrobcu	

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This User Manual has been translated for your convenience using machine translation. Reasonable efforts have been made to provide an accurate translation; however, no automated translation is perfect nor is it intended to replace human translators. The official User Manual is the English version. Any discrepancies or differences created in the translation are not binding and have no legal effect for compliance or enforcement purposes. If any questions arise related to the accuracy of the information contained in the User Manual, please refer to the English version of those contents which is the official version.

Technical data

Parameter description	Parameter value
Product name	Massage chair
Model	PHY_SMC_01
Rated voltage [V~] / frequency [Hz]	230/50
Rated power [W]	90
Max Load [kg]	100
Timer [min]	15-30
Dimensions [width x depth x height; mm]	150 x 108 x 73
Weight [kg]	55.3

1. General description

The user manual is designed to assist in the safe and trouble-free use of the device. The product is designed and manufactured in accordance with strict technical guidelines, using state-of-the-art technologies and components. Additionally, it is produced in compliance with the most stringent quality standards.

DO NOT USE THE DEVICE UNLESS YOU HAVE THOROUGHLY READ AND UNDERSTOOD THIS USER MANUAL.

To increase the product life of the device and to ensure trouble-free operation, use it in accordance with this user manual and regularly perform maintenance tasks. The technical data and specifications in this user manual are up to date. The manufacturer reserves the right to make changes associated with quality improvement. The device is designed to reduce noise emission risks to a minimum, taking into account technological progress and noise reduction opportunities.

Legend

The product satisfies the relevant safety standards.

Read instructions before use.

The product must be recycled.

WARNING! or CAUTION! or REMEMBER! Applicable to the given situation.
(general warning sign)

ATTENTION! Electric shock warning!

Only use indoors.

PLEASE NOTE! Drawings in this manual are for illustration purposes only and in some details may differ from the actual product.

2. Usage safety

ATTENTION! Read all safety warnings and all instructions. Failure to follow the warnings and instructions may result in electric shock, fire and/or serious injury or even death.

The terms "device" or "product" are used in the warnings and instructions to refer to:

Massage chair

2.1. Electrical safety

- a) The plug must fit the socket. Do not modify the plug in any way. Using original plugs and matching sockets reduces the risk of electric shock.
- b) Avoid touching earthed elements such as pipes, heaters, boilers and refrigerators. There is an increased risk of electric shock if the earthed device is exposed to rain, comes into direct contact with a wet surface or is operating in a damp environment. Water getting into the device increases the risk of damage to the device and of electric shock.
- c) Do not touch the device with wet or damp hands.
- d) Use the cable only for its designated use. Never use it to carry the device or to pull the plug out of a socket. Keep the cable away from heat sources, oil, sharp edges or moving parts. Damaged or tangled cables increase the risk of electric shock.
- e) Do not use the device if the power cord is damaged or shows obvious signs of wear. A damaged power cord should be replaced by a qualified electrician or the manufacturer's service centre.
- f) To avoid electric shock, do not immerse the cord, plug or device in water or other liquids. Do not use the device on wet surfaces.
- g) Do not use in very humid environments or in the direct vicinity of water tanks.
- h) Before the first use, please check whether the main voltage type and current comply with the indicated data on the type plate.

2.2. Safety in the workplace

- a) Make sure the workplace is clean and well lit. A messy or poorly lit workplace may lead to accidents. Try to think ahead, observe what is going on and use common sense when working with the device.
- b) Only the manufacturer's service centre may make repairs to the product. Do not attempt to make repairs yourself!
- c) In case of fire, use a powder or carbon dioxide (CO₂) fire extinguisher (one intended for use on live electrical devices) to put it out.
- d) Please keep this manual available for future reference. If this device is passed on to a third party, the manual must be passed on with it.
- e) Keep packaging elements and small assembly parts in a place not available to children.

Remember! When using the device, protect children and other bystanders.

2.3. Personal safety

- a) The device is not designed to be handled by persons (including children) with limited mental and sensory functions or persons lacking relevant experience and/or knowledge unless they are supervised by a person responsible for their safety or they have received instruction on how to operate the device.
- b) To prevent the device from accidentally switching on, make sure the switch is on the OFF position before connecting to a power source.
- c) The device is not a toy. Children must be supervised to ensure that they do not play with the device.

2.4. Safe device use

- a) Do not use the device if the „ON/OFF“ switch does not function properly (does not switch the device on and off). Devices which cannot be switched on and off using the ON/OFF switch are hazardous, should not be operated and must be repaired.
- b) When not in use, store in a safe place, away from children and people not familiar with the device who have not read the user manual. The device may pose a hazard in the hands of inexperienced users.
- c) Keep the device out of the reach of children.
- d) Device repair or maintenance should be carried out by qualified persons, only using original spare parts. This will ensure safe use.
- e) To ensure the operational integrity of the device, do not remove factory-fitted guards and do not loosen any screws.
- f) Do not leave this appliance unattended while it is in use.
- g) Clean the device regularly to prevent stubborn grime from accumulating.
- h) The device is not a toy. Cleaning and maintenance may not be carried out by children without supervision by an adult person.
- i) It is forbidden to interfere with the structure of the device in order to change its parameters or construction.
- j) Keep the device away from sources of fire and heat.
- k) For the following patients, consult your doctor before use:
 - Patients with malignant tumours.
 - Patients with heart attack.
 - Patients with severe dermatosis.
 - Patients with severe bone growths.
 - Pregnancy.
 - Persons undergoing therapy as prescribed by a doctor.
 - Accident-causing abnormalities of the spine or bending of the back.
 - Disease deforming the joints.
 - Persons who have had a pacemaker implanted.
 - Persons with visceral diseases such as gastritis, inflammatory bowel disease, hepatitis, etc.
 - Severe circulatory disorders.
 - People with high fever.
 - Patients with impaired perception.
 - Watch out for children and do not allow them to play on the massage chair.
 - This product is not intended for use by persons (including children) who have limb insufficiency, sensory or mental impairment, or lack the necessary experience and knowledge, unless someone responsible for them.
- l) The recommended operating time is 30 minutes. It is not recommended that the user exceeds the suggested time.

ATTENTION! Despite the safe design of the device and its protective features, and despite the use of additional elements protecting the operator, there is still a slight risk of accident or injury when using the device. Stay alert and use common sense when using the device.


3. Use guidelines

A massage chair is a device that provides massage to the entire body or selected parts of it. Thanks to the various massage programs, the massage chair can be adjusted to your needs and preferences. The massage chair can help alleviate muscle pain, reduce stress, improve blood circulation, as well as relaxation.


The product is intended for home use only!

The user is liable for any damage resulting from unintended use of the device.

3.1. Device description

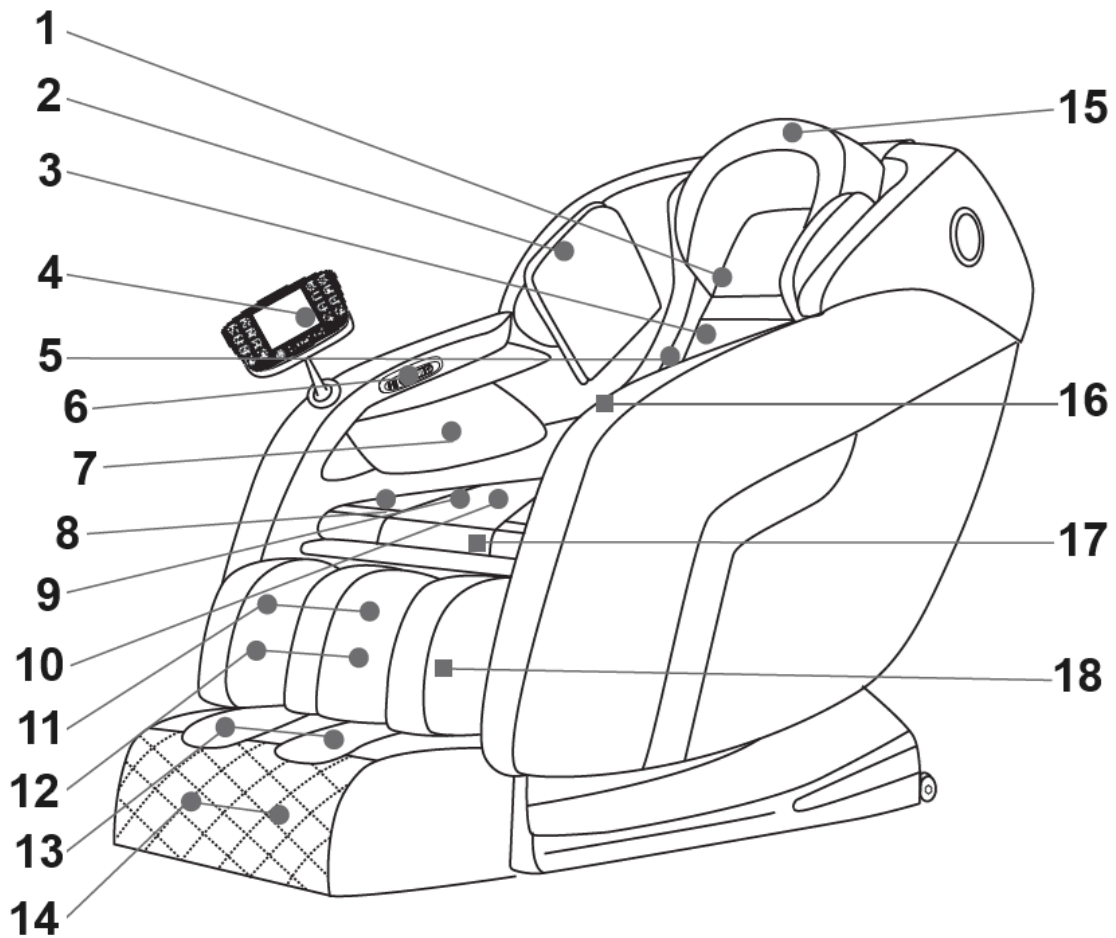1.	Neck Massage	2.	Shoulder Massage
3.	Pillow Pad	4.	Controller
5.	Waist Massage	6.	Quick Controller
7.	Shoulder Airbag Massage	8.	Seat Airbag Massage
9.	Seat Roller Massage	10.	Seat Heating Massage
11.	Calf Airbag Massage	12.	Calf Heating Massage
13.	Calf Heating Massage	14.	Foot Roller Massage
15.	U Shape Pillow	16.	Armrest
17.	Seat	18.	Calf

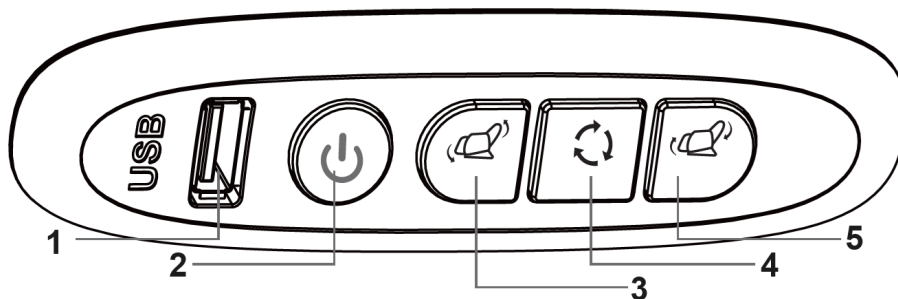- 1 – USB socket
- 2 – On/Off
- 3 – Backrest Up
- 4 – Auto Mode
- 5 – Backrest Down

3.2. Preparing for use


APPLIANCE LOCATION

The temperature of environment must not be higher than 40°C and the relative humidity should be less than 85%. Ensure good ventilation in the room in which the device is being used. There should be at least 40 cm distance between each side of the device and the wall or other objects. The device should always be used when positioned on an even, stable, clean, fireproof and dry surface, and be out of the reach of children and persons with limited mental and sensory functions. Position the device such that you always have access to the power plug. The power cord connected to the appliance must be properly grounded and correspond to the technical details on the product label.

3.3. Assembling the device / Device use

Press the up and down button to position your body in the most comfortable position.

When adjusting the backrest and foot section, make sure there is enough space in the front and back of the seat.

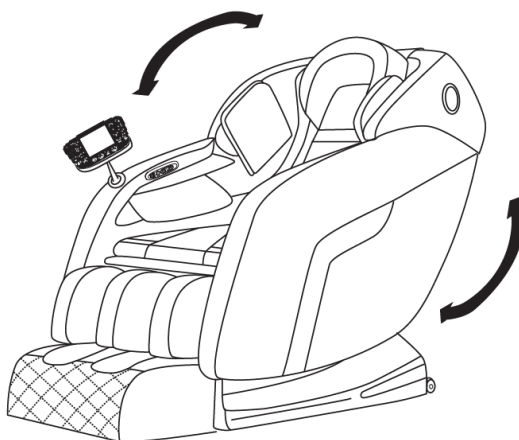3.3.1. How to move the chair

Push the backrest of the massage chair to an angle of approx. 45°.

This will enable the chair to be easily moved to another location.


CAUTION:

- Do not move the massage chair when it is occupied or in use.
- Moving the massage chair may damage the floor. Floor protection is recommended.

3.3.2. How to position the device

Sufficient space must be left at the front and back of the chair.

- Where the massage chair is placed, there must be enough space to recline the seat and raise the foot section.
- Make sure there is at least 40 cm of free space behind the massage chair and at least 40 cm of free space in front of the massage chair.
- Make sure there are no obstructions in front and behind the massage chair.
- Place a rug under the massage chair to prevent damage to the floor.

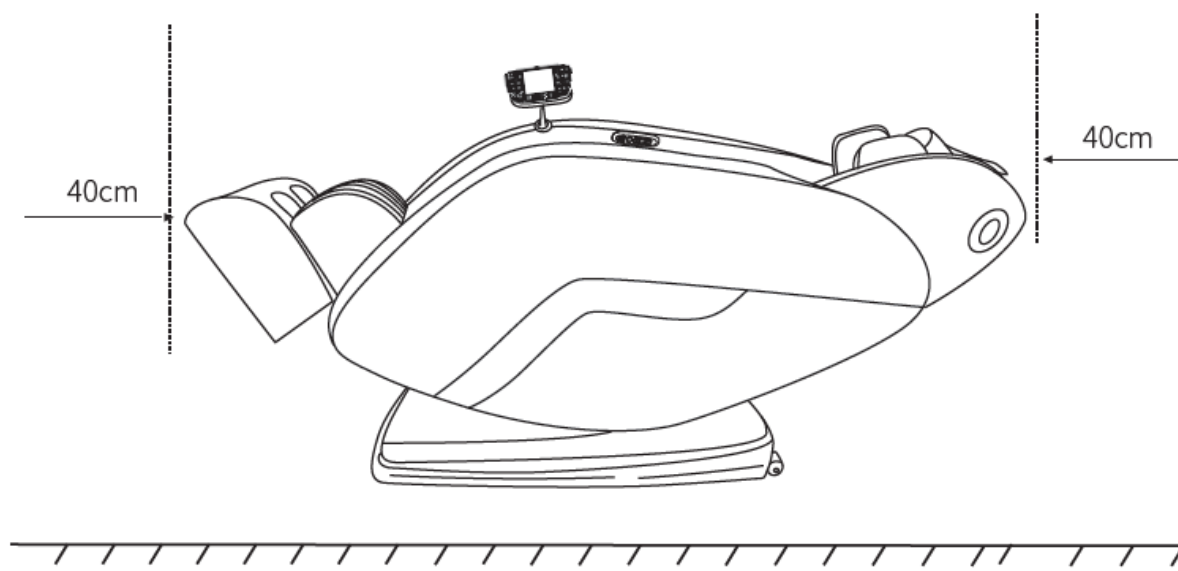
3.3.3. Power ON

Make sure there are no obstructions around the massage chair.

Plug the power plug into the power socket and turn on the power switch.

- Check the power switch before plugging the massage chair into an outlet.
- The power switch should be in the OFF position every time you unplug or plug the chair.
- Follow all instructions on how to properly connect the earth wire.
- Before switching the chair on, fully insert the plug into the socket.

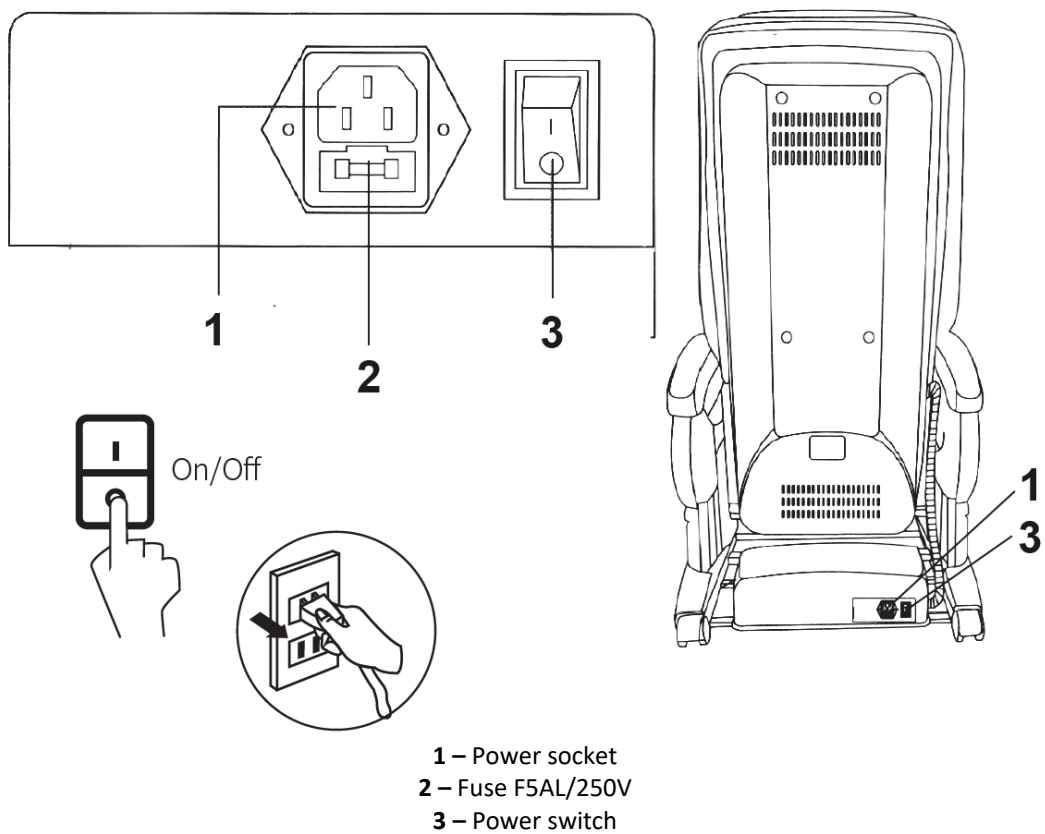1. **Neck** Massage: Press this button to start/stop the neck massage.
2. **Back** Massage: Press this button to start/stop the back massage.
3. **Neck Direction**: Press this button to choose the neck massage direction.
4. **Back Direction**: Press this button to choose the back massage direction.
5. **Head/Shoulder** Airbag Massage: Press this button to start/stop head and shoulder airbag massage.
6. **Arm/Seat** Airbag Massage: Press this button to start/stop arm and seat airbag massage.
7. **Intensity+**: Under manual mode, press this button to increase the massage speed.
8. **Intensity-**: Under manual mode, press this button to decrease the massage speed.
9. **Calf** Airbag Massage: Press this button to start/stop calf airbag massage.
10. **Foot Roller** Massage: Press this button to start/stop the foot roller massage.
11. **Air Intensity+**: Press this button to increase the airbag intensity.
12. **Air Intensity-**: Press this button to decrease the airbag intensity.
13. **Zero Gravity**: Press this button, the chair will go to zero-gravity position.
14. **Roller Speed**: Press this button to switch the speed of the foot roller each time, there are 3 massage speed levels to choose from.
15. **Timing**: Press this button to set the timing of 15-20-25-30 minutes.
16. **Heat**: Press this button to start/stop the heating in the back, seat and calf.
17. **On/Off**: Press this button to start/stop the massage system.
18. **Auto Mode**: Press this button to switch the automatic program each time, there are 20 automatic programs to choose from.
19. **Up**: Press this button, the backrest will be moved up to the original position, release this button, it will stop.
20. **Down**: Press this button, the backrest will be moved down to the lowest position, release this button, it will stop.

3.3.4. Mode Switching

Auto Mode: In auto mode, there are 20 auto programs for choice.

Manual Mode: In manual mode, choose any of manual function to have a specific massage.

3.3.5. Bluetooth connection

The Bluetooth function is activated when the device is switched on.

Open your mobile phone's Bluetooth function to find the device (LEK-988AMY) and connect to it. Open the music player to start playing music.

3.4. Cleaning and maintenance

- a) Always unplug the device before cleaning or putting it away.
- b) Use only non-corrosive cleaners to clean the surface.
- c) After cleaning the device, all parts should be dried completely before using it again.
- d) Store the unit in a dry, cool place, free from moisture and direct exposure to sunlight.
- e) Do not spray the device with a water jet or submerge it in water.
- f) Do not allow water to get inside the device through vents in the housing of the device.
- g) Clean the vents with a brush and compressed air.
- h) The device must be regularly inspected to check its technical efficiency and spot any damage.
- i) Do not leave the battery in the device if it will not be used for a longer period of time.
- j) Use a soft, damp cloth for cleaning.
- k) Do not use sharp and/or metal objects for cleaning (e.g. a wire brush or a metal spatula) because they may damage the surface material of the appliance.
- l) Do not clean the device with an acidic substance, agents of medical purposes, thinners, fuel, oils or other chemical substances because it may damage the device.

DISPOSING OF USED DEVICES:

Do not dispose of this device in municipal waste systems. Hand it over to an electric and electrical device recycling and collection point. Check the symbol on the product, instruction manual and packaging. The plastics used to construct the device can be recycled in accordance with their markings. By choosing to recycle you are making a significant contribution to the protection of our environment. Contact local authorities for information on your local recycling facility.

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